




Absolute Mag™ Anti-CD45 Magnetic Particles Conjugation kit, 200-500 nm

Cat.No: WHM-N038K

DESCRIPTION

Description

Absolute Mag™ Anti-CD45 Magnetic Particles are ideal for rapid isolation or depletion of human CD45+ progenitor stem cells. Absolute Mag™ magnetic particles coated with anti-CD45 monoclonal antibody recognize and efficiently bind to CD45+ cells following a short incubation. The nanoparticle bound CD45+ cells can be separated from the rest of the sample by magnet.

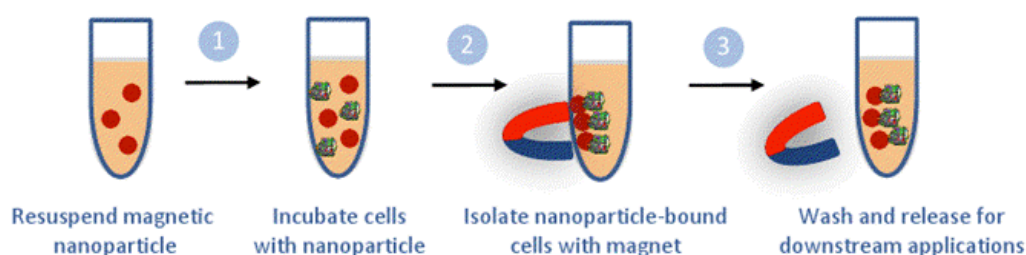

Absolute Mag™ - Anti-CD45 Magnetic Particles offer high recovery of high-purity and viable cells for use in further downstream molecular assays. The beads bound cells can be lysed for further protein or nucleic acid analysis. Absolute Mag™ magnetic particles are much smaller than conventional micro-beads. This feature allows for better accessibility of the particles to the antigenic epitope on cell surfaces. In addition, the surfaces of Absolute Mag™ magnetic particles are uniquely coated to reduce non-specific interactions with non-targeted CD45 negative cells.

KIT COMPONENTS

Kit Components

Absolute Mag™ anti-CD45 magnetic particles are provided in 1 ml of phosphate buffered saline (PBS), pH 7.4.

The kit further includes:

- Washing Buffer
- Elution Buffer

SPECIES REACTIVITY

Species Reactivity

The product has been tested for human CD45+ progenitor stem cells. Others not tested.

STORAGE

Storage Condition

All materials except the magnet should be stored at 4°C. When stored as specified the product is stable for six months.



GENERAL PROTOCOL

Deplete or Positively Isolate CD45+ cells

1. Dilute blood with an equal volume of PBS +2 mM EDTA.
2. Slowly layer the diluted blood over the Ficoll-Hypaque solution in a 50-ml conical centrifuge tube. Use 2 ml of Ficoll-Hypaque (not included, product by Sigma Aldrich) per 1 ml of blood.
3. Centrifuge at 400 g for 40 min at room temperature with no brake.
4. Collect the mononuclear cells located at the interface between plasma and the Ficoll-Hypaque and transfer to 15-ml conical tube.
5. Dilute aspirated mononuclear cells with 4 volume of cold PBS.
6. Centrifuge at 400 g at 4°C for 10 min and discard the supernatant.
7. Re-suspend mononuclear cells from 1ml peripheral blood in 1 ml of Ca²⁺ and Mg²⁺ free PBS buffer supplemented with 0.05% BSA and 5 mM EDTA.
8. Wash 25 µl of Absolute Mag™-CD45 beads twice with 200 µl of PBS Buffer. Re-disperse in 25 µl of PBS buffer.
9. Add pre-washed beads to the mononuclear cell solution with 4-5 million cells and incubate with gentle rotation at 4°C for 1-2 hour.
10. Place the tube by a magnet for 5 min.
For depletion: Transfer supernatant to a new tube for further use and discard the beads.
For positive isolation: While the tube is still in the magnet, carefully remove and discard the supernatant.
11. Repeat the above wash step one more time.
12. Wash beads-bound cells in 200 µl of PBS Buffer and pellet down beads-bound cells as above.
13. Re-suspend isolated cells in PBS or preferred cell medium for further use in downstream applications.